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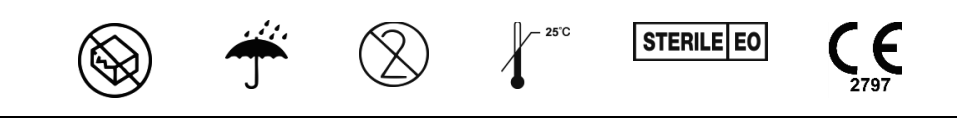
PERFORMANCE. POWER. SPEED.

Device Description

Silverlon® Island Dressings are sterile, multi-layer, absorbent, antimicrobial barrier dressings with attached adhesive tape.

Available as:

Code	Tapе Size	Pad Size
ID23	5 x 8 cm	2.5 x 5 cm
ID34	8 x 10 cm	5 x 6.5 cm
ID44	10 x 10 cm	5 x 5 cm
ID46	10 x 15 cm	5 x 10 cm
ID410	10 x 25 cm	5 x 20 cm
ID412	10 x 30 cm	5 x 25 cm
ID414	10 x 36 cm	5 x 30 cm
ID66	15 x 15 cm	10 x 10 cm

Silverlon® Island Dressings contain silver to provide an ancillary, local antimicrobial effect to reduce the risk of wound infection by organisms sensitive to silver. Ionic silver is a broad spectrum antimicrobial, which has been shown effective against a range of pathogens associated with wound infection. Silverlon® dressings have been tested in vitro and found effective against microorganisms such as: Enterococcus faecalis, Staphylococcus aureus, Pseudomonas aeruginosa, Acinetobacter Baumannii, Enterobacter cloacae, Staphylococcus epidermidis, and Klebsiella pneumoniae.

In vitro testing has shown Silverlon® dressings to be MRI-compatible.

Indications

Silverlon® Island Dressings are professional wound dressings intended for the local of management of:

- Donor and graft site wounds
- Diabetic, pressure and venous ulcers
- Surgical wounds
- Infected wounds*

*Silverlon® Island Dressings may be used on those infected wound types listed above. Where used on infected wounds, the infection should be treated as per local clinical protocol.

Silverlon® Island Dressings are indicated for a total duration of contact (i.e. where treatment involves consecutive application of individual dressings) not to exceed 30 days.

Warnings

- Do not use past expiration date on the product packaging.
- Do not use if pouch is damaged or open.
- Do not use petroleum-based ointments or creams under Silverlon® Island Dressing.
- Do not moisten Silverlon® Island Dressing with hydrogen peroxide or povidone iodine.

Contraindications

- Silverlon® Island Dressings are latex free.
- Avoid using Silverlon® Island Dressings on

Silverlon®

Island Dressings

Silverlon® Island Dressing Instructions for Use

English

- patients with known sensitivity to silver or nylon.
- Silverlon® Island Dressings are not intended for use on 3rd degree burns.

Adverse Reactions

Silverlon® Island Dressings have been subjected to independent standard in vitro and in vivo biocompatibility tests, including cytotoxicity, sensitization and intracutaneous reactivity. All tests were performed in accordance with the International Standard Organization (ISO) 10993 Standard Series for Biological Evaluation of Medical Devices. The results of these studies indicated that Silverlon® Island Dressings are safe for their intended use.

Instructions for Use

- Cleanse wound with sterile water, distilled water, or normal saline, removing necrotic debris or eschar as needed per local protocol.
- Select the Silverlon® Island Dressing with a pad size that overlaps the wound margins by 1-2 cm.
- Activate Silverlon® Island Dressing by thoroughly moistening with sterile water, distilled water, or normal saline; do not moisten the adhesive tape.
- Remove the adhesive tape release liner from one side of the Silverlon® Island Dressing; position the dressing with the pad directly over wound, so that the silver side of the dressing contacts the skin.
- Apply the dressing to the wound and smooth the adhesive tape into place on the intact skin surrounding the wound.
- Remove the release liner from each remaining side of the dressing and smooth the adhesive tape into place on the intact skin surrounding the wound.
- Silverlon® Island Dressings may be used for up to 7 days, but may require more frequent changing depending on wound condition and exudate buildup.
- To remove Silverlon® Island Dressings, gently depress the surrounding skin while lifting the adhesive tape edges.
 - If the dressing pad adheres to the wound, moisten the dressing as needed with sterile water, distilled water, or normal saline, until the dressing can be easily removed.

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